

Updating Luna Glass


Please read the full instructions carefully before starting the update.

To update your Luna Glass, you'll need two things: a USB-C cable and a computer with either Microsoft Edge or Google Chrome installed. **It is currently not possible to update Luna Glass using Safari.**


Connect your Luna Glass to the computer using the USB-C cable, then open the following link:

https://lunaglasses.github.io/LunaWebPageBootloader/user_facing/


The following screen should appear. If any pop-ups requesting permission appear, please click **Allow**. Once this is done, click the blue **Connect** button.

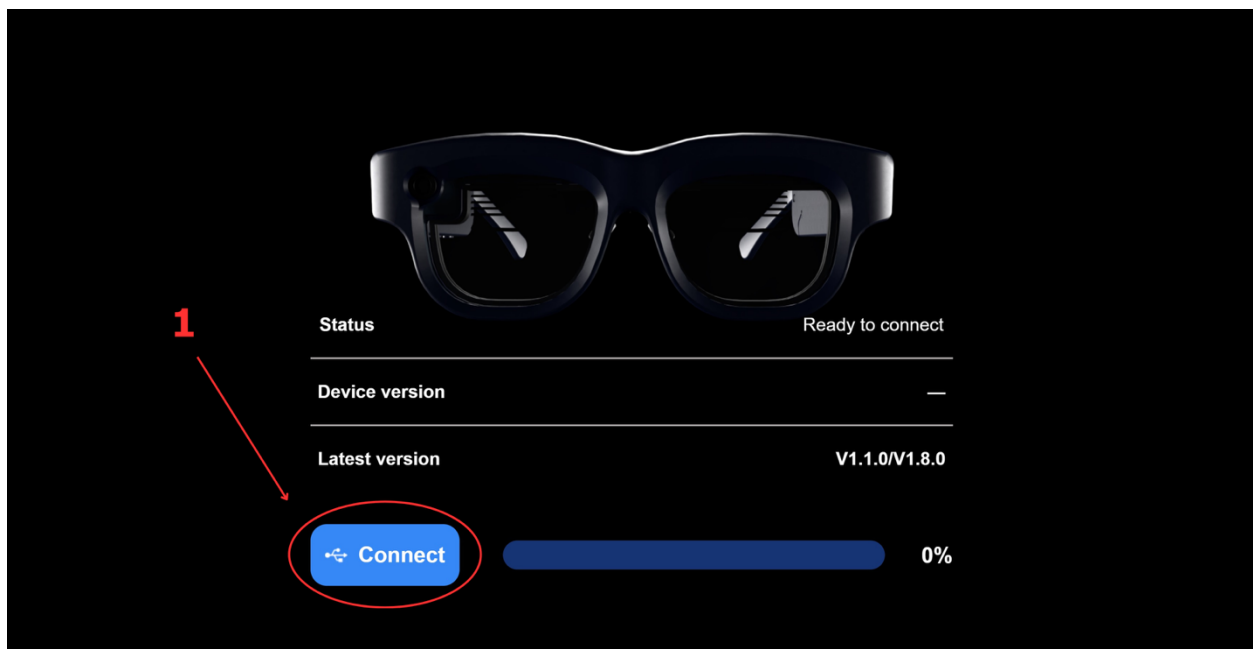

After clicking **Connect**, the following dialog should appear. Select the option that says **Luna Glass Device**, then click **Allow**. This will give the website permission to connect to your Luna Glass.

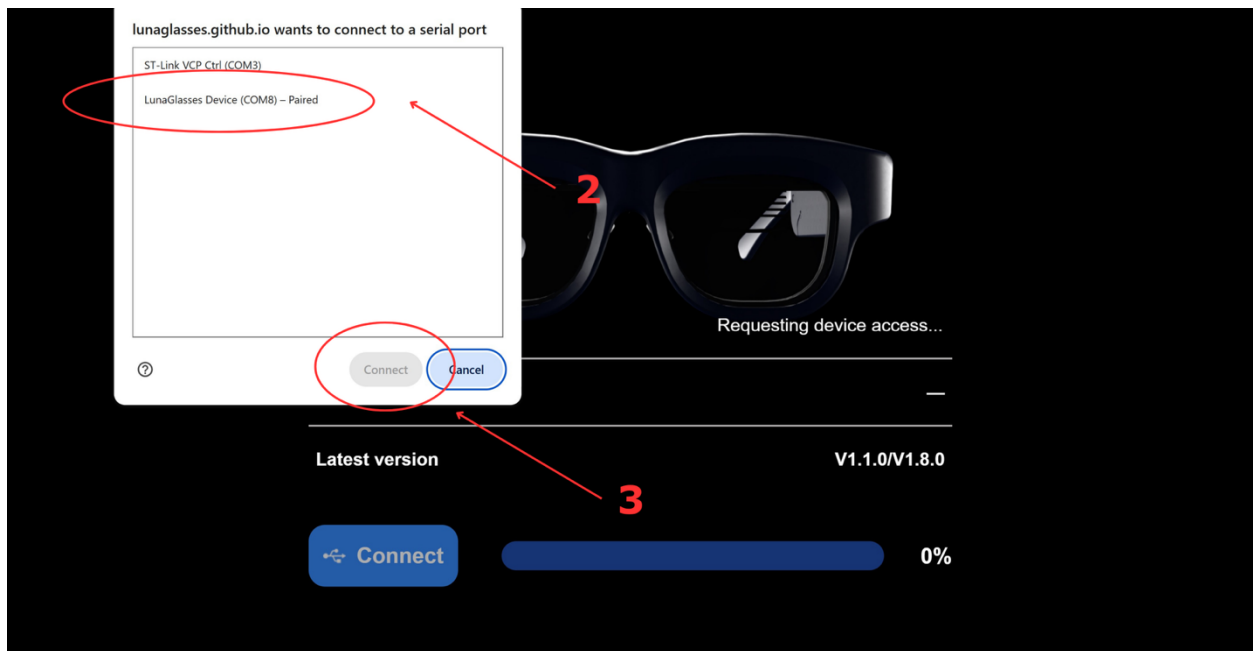


After clicking **Connect**, the website will check whether the Luna Glass firmware is up to date. If an update is available, the option to update the Luna Glass will appear. Click **Update**, as shown below, to start the update process.

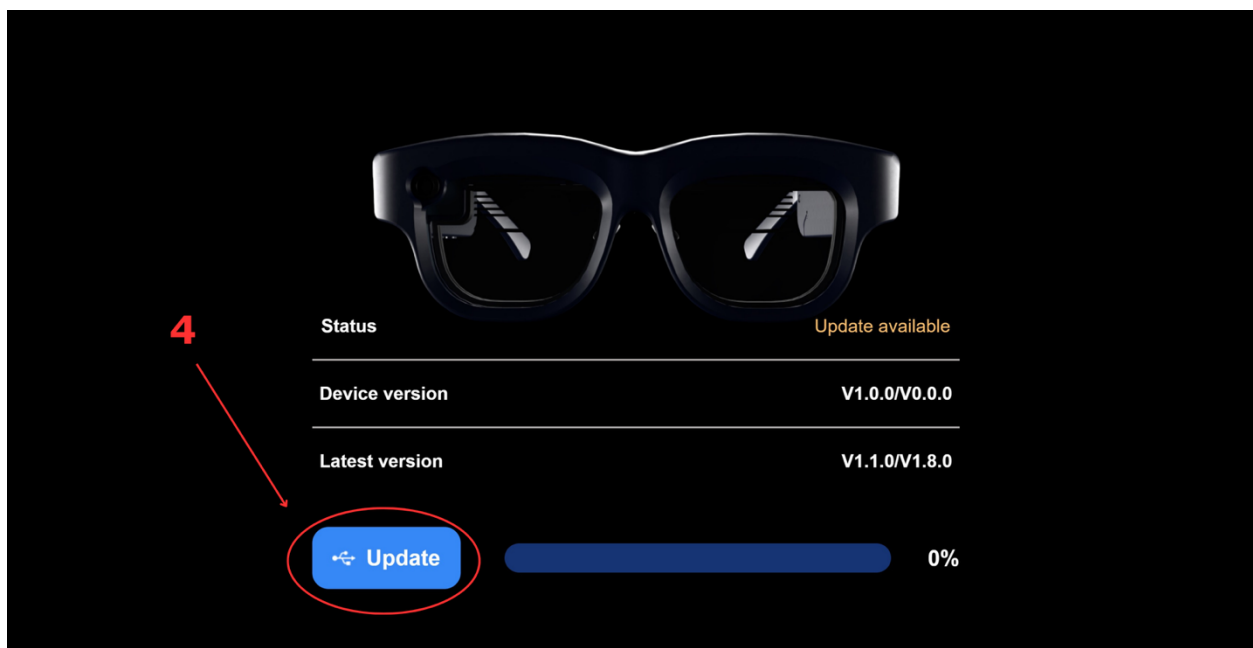

The progress bar will show the status of the update. It is important **not to unplug or turn off the Luna Glass** until the update is complete.




Once the update is complete, you can unplug and restart the Luna Glass. The update is now complete, and your Luna Glass is ready to use!

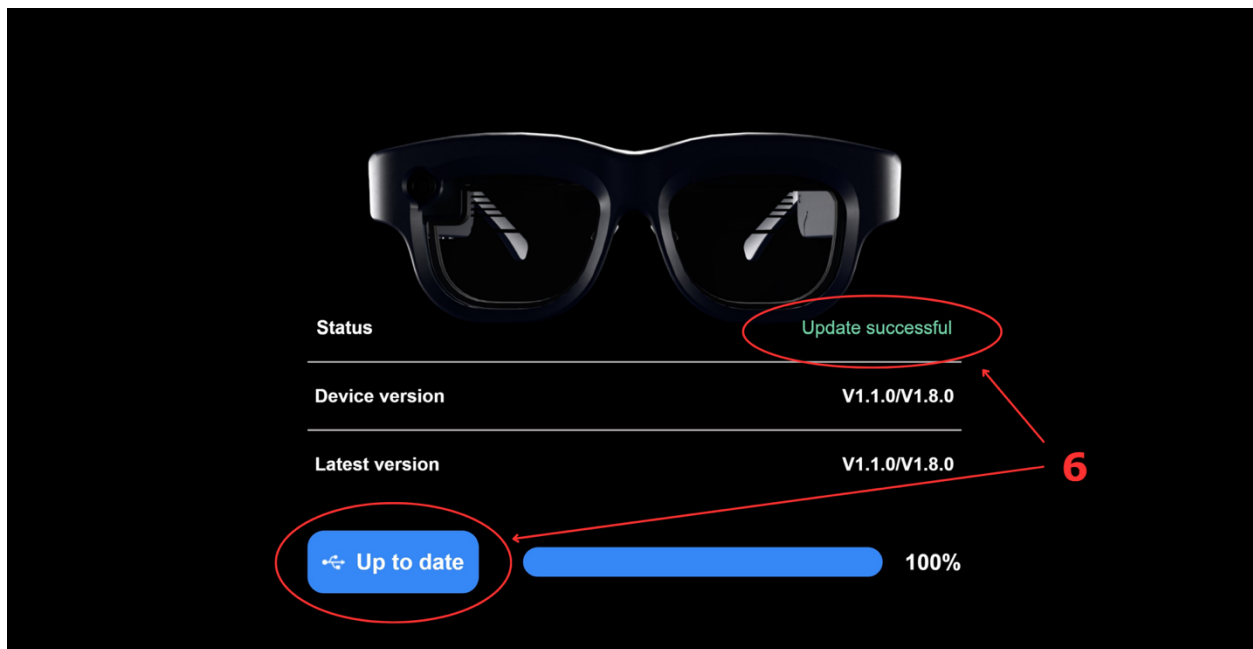

Unauthorised Devices

If the following dialog box appears after connecting the Luna Glass, please copy the **Device ID** and email it to support@lunaglasses.nl.

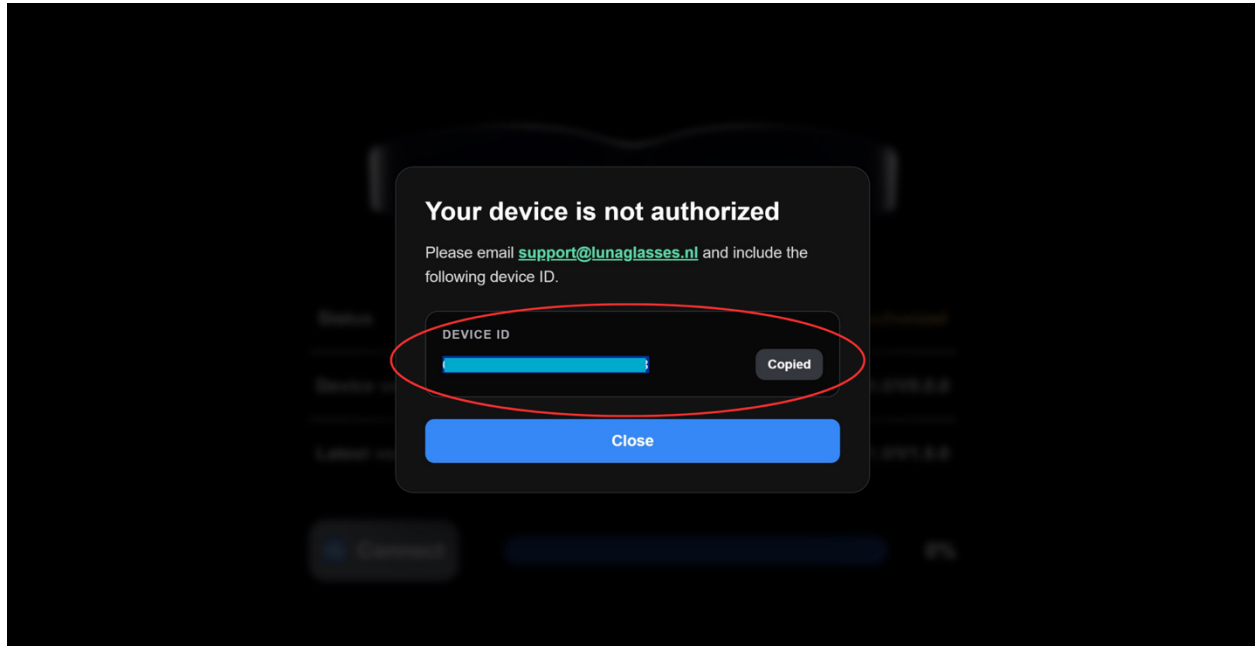

If you don't have **Microsoft Edge** or **Google Chrome** installed, you can download them using the links below:

- [Download Microsoft Edge](#)
- [Download Google Chrome](#)

If you experience any issues during the update process, please contact us at support@lunaglasses.nl